

Ion exchangers for demineralising water.



**Safe, reliable, top quality
that can be regenerated
time and time again.**

stakpure

H₂O – WASSER IN SEINER REINSTEN FORM



Ion exchangers from stakpure. The economic answer to demineralised water.

Whether for ultrapure water systems, automatic analysers, dishwashers, autoclaves – demineralised water is needed on a daily basis for all sorts of different purposes in medical practices, hospitals and laboratories. Whether in small quantities or up to 3000 litres an hour – we have the right ion exchanger for any need.

Inside a container, water flows over a mixed-bed of special resins that removes its salts and other minerals. A tried and proven method that guarantees water of a constant high quality all the time. The resins come with a long useful life and can be regenerated whenever necessary without any impact on the environment.

For small quantities: Ion exchangers in cartridges of synthetic material.

These ion exchangers have been designed specifically for medical practices and small laboratories requiring up to 10 litres a day. The non-pressurised cartridge will be simply regenerated or replaced after a prolonged period of service. It is equipped with conductivity meter, set of hoses and wall mounting is optionally available.

Operating safety

An innovative filter system prevents any resin from escaping. High product safety for downstream equipment as well.

stakpure ion exchangers	DS 450*	DS 750*
Model	one-way cartridge	regenerably cartridge
Flow rate [l/h]	50	100
Capacity for 10° total dissolved solids [Liters]	425	750
Pure water quality [$\mu\text{S}/\text{cm}$]	0.1–20	0.1–20
Max. water temperature [°C]	30	30
Power supply [V/Hz]	230/50–60	230/50–60
Material	Polyethylene	Polyethylene
Dimensions [\varnothing x height mm]	100 x 600	175 x 470
Connection	R3/4"	R3/4"
Weight, approx. [kg]	3	6
Order no.	11042500	11170500
Replacement cartridge	11042700**	11175000
Wall mount	– included	13140400

* conductivity meter included

** set 2 pcs.



For large quantities: Ion exchangers in stainless steel cartridges.

For any situation requiring 300 to 3000* litres of demineralised water an hour.
A sophisticated flow system for incoming water ensures a fast and effective flow
of water through the mixed-bed system with ion-exchange resins.

* We recommend our reverse-osmosis systems in applications requiring higher quantities or
where particular demands are placed on pure water.

stakpure DS ion exchange cartridges


For up to 300 and 1000 l an hour.
Stainless-steel cartridge. Hard-rubber carrying handles.



stakpure mixed-bed ion exchangers	DS 1500	DS 2000	DS 2800	DS 4000	DS 6000
Water demineraliser			pressure-resistant		
Stainless-steel cartridge			pressure-resistant to 10 bar		
Flow rate [l/h]	300	300	950	1000	1000
Capacity for 10° total dissolved solids [Liters]	1500	2000	2800	4000	6000
Pure water quality [µS/cm]			0.1 – 20		
Max. water temperature [°C]			30		
Material			Stainless steel – AISI 316		
Dimensions [Ø x height mm]	240 x 410	240 x 490	240 x 570	240 x 700	240 x 1155
Connection			R3/4"		
Weight, approx. [kg]	14	18	24	27	45
Order no.	12115000	12200000	12280000	12400000	12600000
Wall mount			13140500		

On request, the cartridges can also be equipped with quick-connect couplings.

**Comparison
Pure water quality
in relation to
capacity**

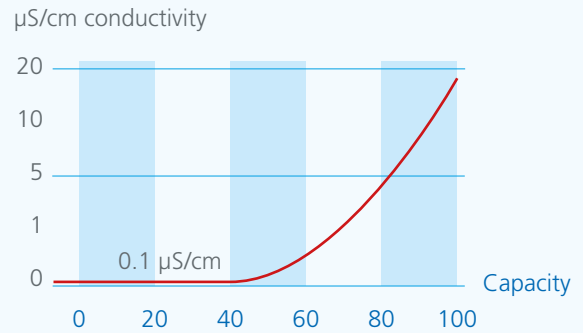

stakpure ion exchangers deliver reproducible demineralised water with a permanently low level of conductivity that remains more or less constant throughout cartridge life and meets all common standards, such as ASTM, CLSI and DIN EN 285/13060.



stakpure industrial ion exchange cartridges

Specifically for industrial applications. With a particularly high capacity and flow rates of up to 3000 litres an hour. Made entirely of stainless steel.

stakpure mixed-bed ion exchangers	DS 7000	DS 11000	DS 15000
Water demineraliser		pressure-resistant	
Stainless-steel cartridge		pressure-resistant to 10 bar	
Flow rate [l/h]	2000	2500	3000
Capacity for 10° total dissolved solids [Liters]	7000	11000	15000
Pure water quality [µS/cm]		0.1–20	
Max. water temperature [°C]		30	
Material		Stainless steel – AISI 316	
Dimensions [Ø x height mm]	363 x 600	363 x 850	363 x 1100
Connection		R3/4"	
Weight, approx. [kg]	55	70	90
Order no.	12700000	12110000	12150000

On request, the cartridges can also be equipped with quick-connect couplings.

Everything under control. stakpure digital and analogue conductivity monitoring.

stakpure DC 400

Developed for online measurement of pure and ultrapure water downstream of ion-exchange cartridges and for monitoring ring-main systems. High-precision conductivity meter ensuring exact measurement of temperature and conductivity for constantly monitoring quality. Digital display guarantees to show readings even at the lowest end of the conductivity measurement range.

Comes complete with connecting cable, plug-type power supply unit, set of two 1.5 m hoses and measuring cell.

Voltage selection

Automatic voltage selection for international use.

LED indicators

Show red/green to signalise when limit values are exceeded for temperature and conductivity.

Technical specifications	DC 400
Conductivity measurement range [$\mu\text{S}/\text{cm}$]	0.1 – 300 (automatically switches over at 9.99)
Temperature measurement range [$^{\circ}\text{C}$]	0.1 – 50
Temperature compensation	automatic
Conductivity limiting-value setting [$\mu\text{S}/\text{cm}$]	0.1 – 30, infinitely variable
Temperature limiting-value setting [$^{\circ}\text{C}$]	1 – 50, infinitely variable
Signals on exceeding limiting value	red/green LED indicator + buzzer (both limiting values can be deactivated)
Fault signalling	Display changes from green to red
Interface (optional only)	potential-free contact
Power supply [V/Hz]	110/230/50–60 (automatic voltage selection)
Dimensions [W x D x H mm]	73 x 32 x 141
Installation	Wall mounting
Connection	R3/4"

Unit	DC 400	DC 400	DC 400
Version		with potential-free contact	with potential-free contact + solenoid valve
Order no.	14180500	14180600	14180700

On request, the conductivity meters can also be equipped with quick-connect couplings.

Large graphic display

2-colour, back-lit graphic digital display indicating conductivity in $\mu\text{S}/\text{cm}$ or in $\text{M}\Omega\text{cm}$ and temperature in $^{\circ}\text{C}$. Background illumination changes from green to red to signal a fault and show plaintext message.



Interfaces

Direct connection of peripheral equipment at potential-free contact.

Limiting values

Infinite adjustment of limiting values for temperature and conductivity.

High-precision measurement

Automatic calibration before each measurement with integrated reference resistor, temperature measured by platinum temperature sensor.

stakpure AC 100

Conductivity meter as the easy and reliable answer to measuring the conductivity of pure water. Analogue display of measurement value in $\mu\text{S}/\text{cm}$. Either for connecting directly to ion-exchange cartridges or for space-saving wall installation.

Comes complete with measuring cell and set of two 1.5 m hoses.

Unit	AC 100	AW 100
Measurement range [$\mu\text{S}/\text{cm}$]	0.1 – 50	0.1 – 50
Power supply [V/Hz]	230/50 – 60	230/50 – 60
Dimensions [W x D x H mm]	93 x 123 x 170	93 x 125 x 96
Connection	R3/4"	R3/4"
Order no. (for connecting directly to DS 1500–15000)	14160100	
Order no. (for wall installation only)		14160200

On request, the conductivity meters can also be equipped with quick-connect couplings.



AC 100



1.
Pure water
outlet valve



2.
Pure water
outlet valve



3.
Pure water
tap



4.
Pure water
distributor block



5.
Pure water
T-connector

System accessories	Connections	Order no.
1. Pure water outlet valve Synthetic, straight	R3/4" with union nut	13140000
2. Pure water outlet valve Synthetic, angled	R3/4" with union nut	13140100
3. Pure water tap Synthetic	R3/4" with union nut	13140300
4. Pure water distributor block Synthetic	R3/4" 1 inlet and 3 outlets for connecting tap and distributing valves	13140200
5. Pure water T-connector Synthetic	R3/4" 1 inlet and 2 outlets for connecting tap and distributing valves	13142000



stakpure GmbH

Auf dem Kesseling 11

56414 Niederahr

Germany

Phone: +49 (0) 2602 10673-0

Fax: +49 (0) 2602 10673-200

info@stakpure.com

www.stakpure.de



We are certified
according to
ISO 9001: 2015

Is reliable and economic preparation of pure and/or ultrapure water a topic for you?
Just call us!

info@stakpure.com
www.stakpure.de

Retailer panel